


# HAND HELD PULSE OXIMETER




## H100B Pulse Oximeter

- LCD display
- RTC (Real Time Clock) display
- Backlight control and automatic power-off function for power saving
- SpO<sub>2</sub> and Pulse rate measurement
- Display numeric and waveform of SpO<sub>2</sub> simultaneously
- Long battery life-up to 48 hours on 4 "AA" or 36 hours on Ni-MH rechargeable batteries
- Powerful data storage capacity
- Up to 300 hours SpO<sub>2</sub> data (1 read per 2 seconds) can be transferred to PC by Patient Care Viewer for storage, review and printing
- The latest 10 minutes trend graph and table of SpO<sub>2</sub> and Pulse Rate can be reviewed in the screen
- Audible and visual alarm
- Pitch tone

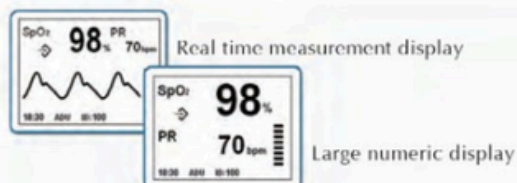


PC storage is easy to access through RTC and ID

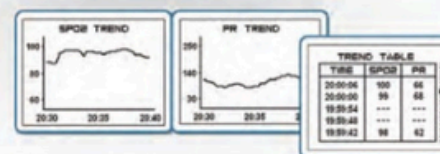

The latest 10 minutes trend graph and table of SpO<sub>2</sub> and Pulse Rate can be reviewed in the screen



Battery charger stand

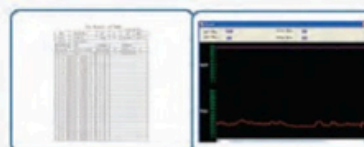

H100B provides up to 100 IDs and 300 hours data storage capacity for continuous monitoring, all data can be transferred to PC through DB9 interface for storage, review and printing.



H100B can support reusable and disposable DB9 sensors for adult, pediatric and neonatal monitoring.